

3 IN 1 PADDED COMMmode

LVA3113

vive
vivehealth.com




To see all FAQs
in one place,
visit [vhealth.link/cu2](https://vivehealth.link/cu2)

OWNER'S MANUAL

3-IN-1 PADDED COMMODE


WHAT'S INCLUDED

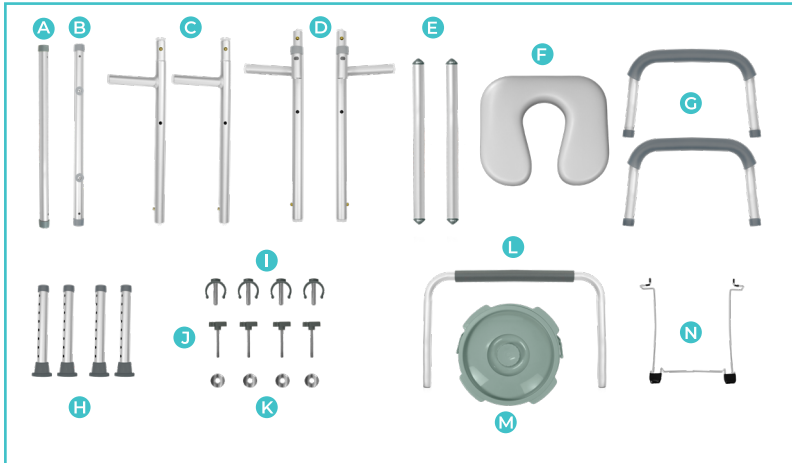


- 1x Front Tube (A)
- 1x Back Tube (B)
- 2x Front Support rail (C)
- 2x Back Support rail (D)
- 2x Reinforced Crossbar (E)
- 1x Padded Seat (F)
- 2x Armrest (G)
- 4x Feet (H)
- 4x C Bolt (I)
- 4x Easy Twist Bolt (J)
- 4x Washer (K)
- 1x Backrest (L)
- 1x Bucket (M)
- 1x Wire Rack (N)

ASSEMBLY INSTRUCTIONS


1. Insert the front support rails (C) into the front tube (A). Secure with C bolts (I).

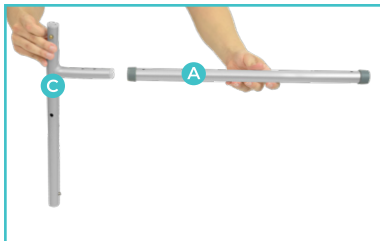


2. Insert the back support rails (D) into the back tube (B). Secure with C bolts (I).

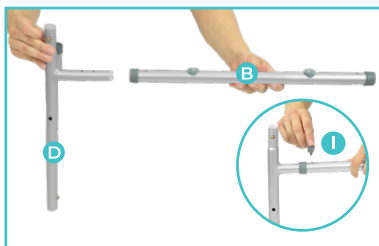

3. Attach all four feet (H) to the front (C) and back support rails (D).

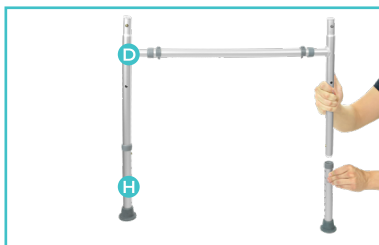

4. Connect the front (C) and back (D) support rails with reinforced crossbars (E). Secure using easy twist bolts (J) and washers (K). Ensure the arrow icon faces upward.

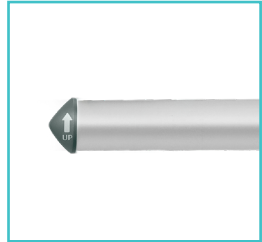

5. Attach both armrests (G) to the frame, ensuring they angle outward.



6. Install the backrest (L) by pressing the push pins and sliding into place.



7. Mount the wire rack (N) into the designated frame openings.



8. Place the padded seat (F) onto the frame.



9. Insert the bucket (M) into the wire rack (N).



WARNINGS ⚠

- Ensure feet are dry and free of soap residue before each use.
- Maximum weight capacity: 330 lbs (150 kg).
- Do not stand on the chair.
- Do not lean excessively on the backrest — this may cause tipping.
- Never attempt to stand from the chair onto a wet or slippery surface.
- Before each use, verify all four legs are locked at equal height.
- Periodically check all screws to ensure they are fully tightened.
- Before each use, verify all four feet are locked at equal height.
- Periodically check all screws to ensure they are fully tightened.

CARE INSTRUCTIONS

- **Daily Cleaning:** After each use, remove the bucket (M), empty contents, and wash thoroughly with mild soap and warm water. Rinse and dry before reinserting.
- **Disinfection:** Regularly disinfect the bucket (M), seat (F), and armrests (G) with a household disinfectant or a solution of 1 part vinegar to 10 parts water.
- **Frame Maintenance:** Wipe the frame with a damp cloth and mild detergent. Avoid abrasive cleaners that may damage the surface.

- **Drying:** Ensure all parts are completely dry before reassembly to prevent rust or mold growth.
- **Inspection:** Periodically check bolts (I, J), washers (K), and feet (H) to ensure they are secure and free of wear.
- **Storage:** Store the commode in a clean, dry area when not in use. DO Not expose to prolonged direct sunlight.

WARRANTY INFORMATION

Vive Health Limited Warranty – 60 Days

Vive Health offers a 60-day limited warranty on the 3-in-1 Padded Commode. This warranty covers defects in materials and workmanship under normal use and conditions.

What is Covered:

- Manufacturing defects or material failures.
- Replacement or repair at no cost if the product is deemed defective.

What is Not Covered:

- Normal wear and tear.
- Damage resulting from misuse, improper care, or unauthorized modifications.
- Accidental damage such as cracks, dents, or stains.
- Use outside of the product's intended function.

Submitting a Warranty Claim:

To initiate a claim, contact Vive Health Customer Support at:

- Web: service@vivehealth.com
- Tel. no. 1-800-487-3808

Provide the following:

- Proof of purchase (e.g., order confirmation or receipt)
- A clear description of the issue
- Photos or videos supporting your claim

NOTES



HAVE MORE QUESTIONS?

Check out our list of Frequently Asked Questions at
vhealth.link/cu2 for helpful answers.



And if that doesn't answer your question, our
customer service team would love to help!
Feel free to connect with them by phone,
e-mail, or chat on our website

service@vivehealth.com
1-800-487-3808
vivehealth.com

Distributed by



8955 Fontana Del Sol Way
Naples, FL 34109
Made in China